
Caloric Irrigators Air or Water



The irrigators from Micromedical by Interacoustics offer a choice of stimuli, either water or air, for caloric testing during the VNG test protocol




Micromedical
by Interacoustics

Air Fx Caloric irrigator

Highlights:

- Cooling and heating ability: 20°C to 50°C
- Illuminated and magnified otoscope
- Self-contained cooling ability
- Integrated with our VNG software

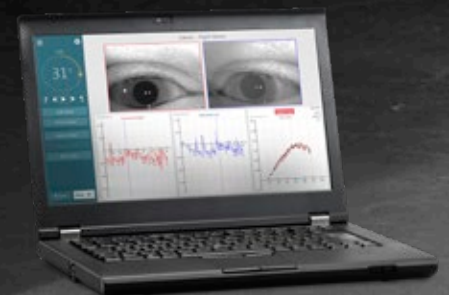

Caloric test with the VisualEyes™ 515/525 software

Versatile Caloric Air Irrigator

The AirFx caloric irrigator has many benefits. It offers an otoscopic handle with illumination as well as magnification. In addition, irrigations can be started using the start button on the handle or via the VisualEyes™ 515/525 or VN415/VO425 software.


The Air Fx is available in 110-130 V and 220-240 V configurations, and is available world-wide.

An integrated water reservoir allows for cool air irrigations down to 20°C regardless of ambient temperature, making it suitable for warm climate use. For warm air irrigations, the Air Fx is capable of heating air up to 50°C.

Temperature Accuracy and Stability is within +/- 0,5°C, and air flow stability is within +/- 0,5 L/min.



Aqua Stim Caloric irrigator

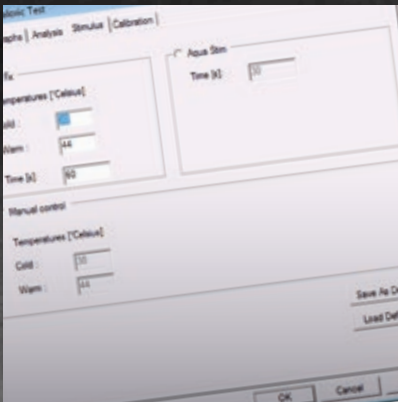

Can be operated from the VisualEyes™ 515/525 or VN415/VO425 software.

Versatile Caloric Water Irrigator

The Interacoustics Aqua Stim offers reliable caloric irrigation. The Aqua Stim can be used as a stand-alone device, or you can control it from the VisualEyes™ 515/525 or VN415/VO425 software.

Irrigation is started by simply pressing the start button on the handle one time to initiate the user-defined caloric protocol.

An external water tank has capacity for approximately 14 irrigations, and eliminates the need for a drain in the examination room.

The Aqua Stim is available world-wide and is supplied as a 110-130 V or 220-240 V model.


Highlights:

- 30°C cold and 44°C warm irrigations
- Use with our VNG software or as stand-alone
- LED illumination and start button in handle
- Eliminates the need for drain and faucet in examination room



Included external water tank



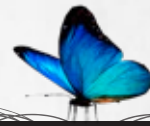

Read more at
micromedical.com
interacoustics.com

academy.interacoustics.com (training materials)

Go online to
explore our
full product
range



Micromedical

mminfo@micromedical.com
800-334-4154
217-483-2122
micromedical.com

Micromedical
by Interacoustics